



WESTERN BLOT > Owi ECL RevelBIOT® Femto

Reliability. Accuracy. Success.


We are with you to improve your research.

PCR

Polymerase (Standard, Hot Start, High Fidelity), dNTP mix

RT-qPCR

Master Mix, dNTP mix

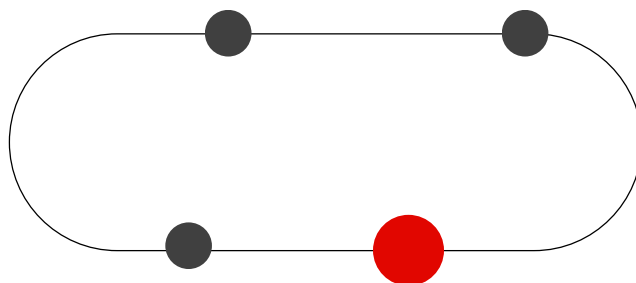

ELECTROPHORESIS


*Agarose, Dye DNA Stain,
DNA Ladder*

WESTERN BLOT

- **Owi ECL RevelBIOT® Femto 250 ml**
- Owi ECL RevelBIOT® Pico 100 ml
- Owi ECOlor® Prestained Protein Marker

Owi ECL RevelBIot® Femto


- **High sensitivity** : Mid to low femtogram detection
- **High intensity** combined with **long-lasting chemiluminescence**
- **Working solution stable** for at least three days

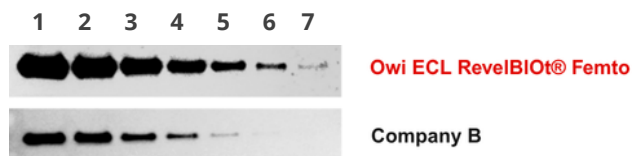

Serial dilution of the protein fraction of a cell extract:

1: 5,000 μ g
2: 2,500 μ g
3: 1,250 μ g
4: 0,625 μ g
5: 0,313 μ g
6: 0,156 μ g
7: 0,078 μ g

Extended signal duration with Owi ECL RevelBIot® Femto

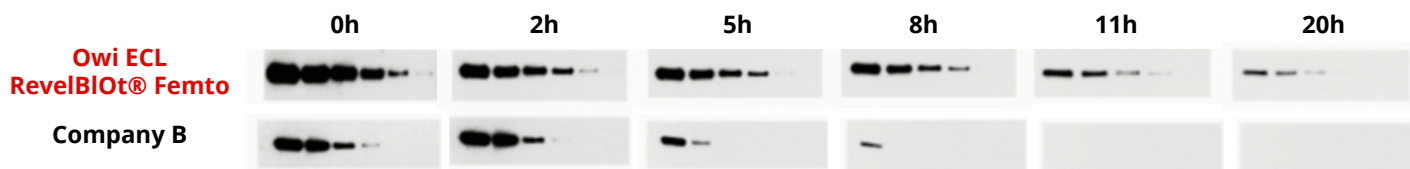

Reference	Name	Size
OZYB002-1000	Owi ECL RevelBIot® Femto	100 ml

[VISIT OUR WEBSITE TO ORDER](#)

Pick the Right ECL RevelBIOT®

	Owi ECL RevelBIOT® Pico	Owi ECL RevelBIOT® Femto
Kit features	Routine Analysis	Proteins with low to medium expression
Sensitivity	> 1 - 10 pg	< 1 pg
Primary antibody	1:500 - 1:5000 0,2 - 2 µg/ml	1:1000 - 1:50 000 20 ng -1 µg/ml
Secondary antibodies	1:20 000 - 1:100 000 10 to 50 ng/ml	1:50 000 - 1:250 000 4 to 20 ng/ml
Exposure time	< 2 min	> 5 min
Other ECL substrates	Pierce ECL, GE ECL, Pierce SuperSignal Pico, GE ECL Plus	Pierce SuperSignal Dura, GE ECL Advance, Pierce SuperSignal Femto
Signal duration	8 hours	24 hours

[VISIT OUR WEBSITE TO ORDER](#)

Your local vendor :